

GEI Consultants Golden Algae Screen

Prepared: October 6, 2025

Prepared By: GreenWater Laboratories/Andrew Chapman

<u>Sample ID</u>	<u>Site</u>	<u>Collected</u>
CL-10	Canyon Lake	30 September 2025

Method

One mL aliquot of live CL-10 sample was transferred to an Utermöhl Counting Chamber. Sample was scanned for the presence of Golden Algae (*Prymnesium parvum*) using a Nikon Eclipse Ti-2 Inverted Microscope equipped with phase contrast optics and epifluorescence. A second one mL aliquot was lightly preserved with Lugol's Iodine solution, allowed to settle and 50 random fields were observed at 600X to estimate *P. parvum* densities.

CL-10

Estimated cell density of *Prymnesium parvum* (Fig. 1) was 1,005 cells/mL.

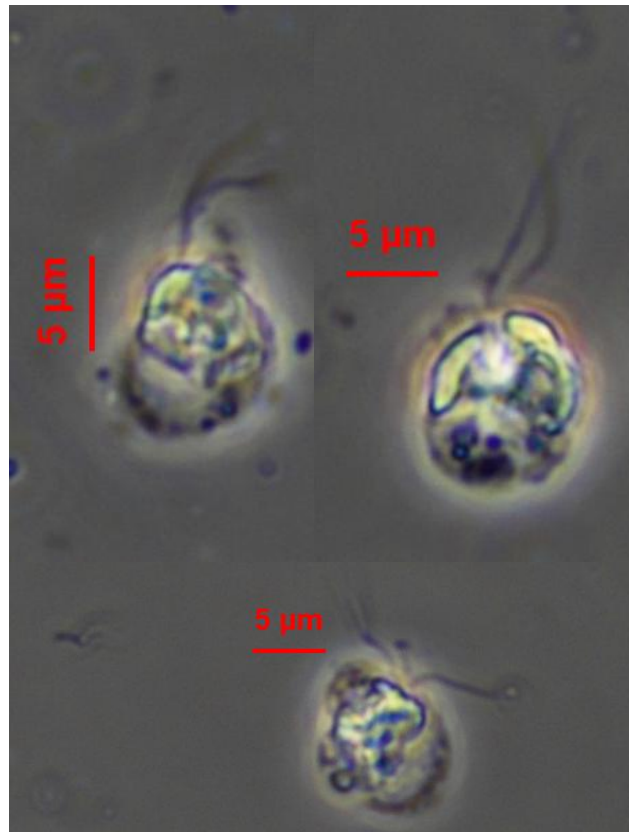


Fig. 1 *Prymnesium parvum* 600X Lugol's Preserved

Submitted by:



Date:

October 6, 2025

Andrew D. Chapman

The results in this report relate only to the samples listed above.

This report shall not be reproduced except in full without written approval of the laboratory.